

MAHURKAR* Acute Triple Lumen Catheter

THE TRIPLE LUMEN
CATHETER
INDICATED FOR:

HEMODIALYSIS

APHERESIS

INFUSION




ARTERIAL LUMEN SLOT DESIGN

Helps reduce positional occlusion while indwelling.



MODIFIED 'DOUBLE D*' LUMEN DESIGN

Ensures high flows with low pressure.

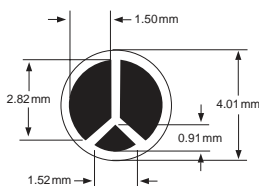

CONVENIENT CATHETER KITS

MAHURKAR* Insertion Kits simplify stocking and inventory tracking and provide most items needed for catheter insertion.

CATHETER KIT COMPONENTS

- (1) 12Fr. catheter
- (1) 18 gauge introducer needle
- (1) J/Straight 0.035" guidewire
- (1) 10Fr. dilator
- (1) 12Fr. dilator
- (2) adhesive wound dressings
- (3) sealing caps



SELECTION

- The 12Fr. MAHURKAR* Acute Triple Lumen Catheter is available in four lengths: 13cm, 16cm, 20cm and 24cm to accommodate varying body types and insertion sites
- Catheters are available individually or in kits

SAFETY

- Soft Atraumatic Tip designed to reduce vessel trauma during insertion
- Laser Cut Arterial Slot reduces positional occlusion, maintains high flows and retains heparin lock
- Imprinted Priming Volumes on catheter extensions quickly identify heparin loading volumes
- Straight Infusion Lumen differentiates it from the dialysis ports to prevent accidental infusion of heparin or Urokinase

EASE OF USE

- 180° Curved Extensions help stabilize the catheter and aid in patient comfort
- Rotatable Suture Wings enable the tip to be rotated after implantation and stabilize the catheter

Ordering Information

REF	Product Description	Units per Box
Catheter Kits		
8888345603	12Fr. x 13cm	5
8888345611	12Fr. x 16cm	5
8888345629	12Fr. x 20cm	5
8888345637	12Fr. x 24cm	5
Catheters		
8888340629	12Fr. x 20cm	5
8888340637	12Fr. x 24cm	5
8888345512	12Fr. x 16cm	5
8888345520	12Fr. x 20cm	5
8888345538	12Fr. x 24cm	5

* MAHURKAR and DOUBLE D are a trademark of Sakharam D. Mahurkar, MD and are used under license.
© 2003 Tyco Healthcare Group LP. All rights reserved.